




# Applicator


For application of Allflex A-Tag™ one-piece ear tags  
A-Tag™ Application Instructions




Flip applicator pin out to load A-Tag™

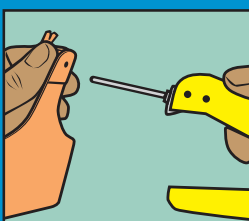


Fully insert applicator pin into A-Tag™ while rotating tag anchor

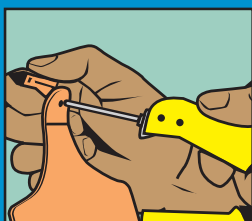


Seat tag anchor "v" onto pin so that pointed tip flexes outward

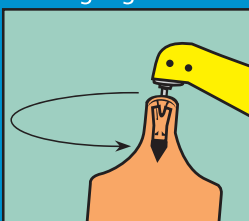

Rotate pin and A-Tag™ back into tagging position

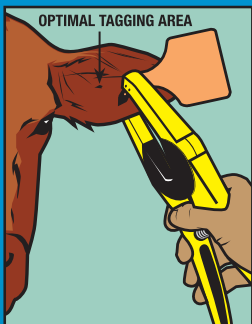

Position tag point between the two cartilage ribs in ear. Do not use existing ear tag holes

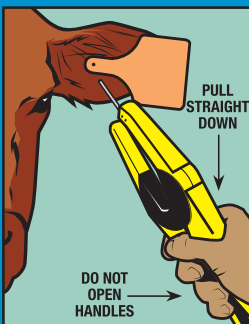

To apply tag, squeeze applicator handles together. Pull applicator straight down, away from animal



ALLFLEX USA, INC.  
P.O. Box 612266  
Dallas/Fort Worth Airport, TX 75261-2266  
(800) 989-TAGS (8247)  
[www.allflexusa.com](http://www.allflexusa.com)

ALLFLEX CANADA  
4135 Bérard  
St-Hyacinthe, QC J2S 8Z8  
866-505-TAGS